
medisana®

BS 483

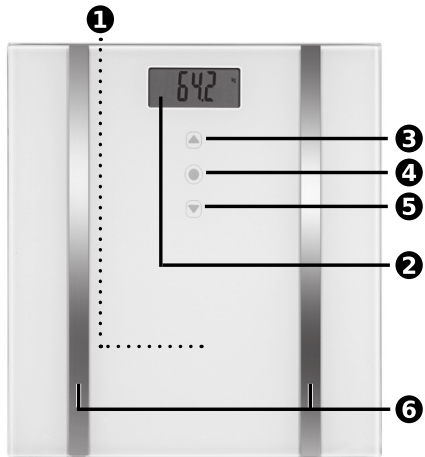medisana
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

REF 40449

EAN 4015588 40449 8

Körperanalysewaage BS 483

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Body Analysis Scale BS 483

- ① Battery compartment (on underside)
- ② Display
- ③ ▲ Up button
- ④ SET button
- ⑤ ▼ Down button
- ⑥ Electrodes

Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- These scales are not suitable for people with pacemakers or other medical implants. The body fat information may be inaccurate in people with diabetes or other medical/physical limitations. The same also applies to very highly trained athletes.
- This device is not suitable for pregnant women!
- The scale has been manufactured for domestic use. It is not suitable for commercial use in hospitals or other medical institutions.
- Do not use the scale in a location where extremely

high temperatures or humidity may occur. Keep the scale away from water.

- Do not place the scale on an uneven surface. Do not stand on the corners/edges of the scale. When weighing yourself, do not stand on one side or the edge of the scales.
- Danger of slipping! Do not stand on the scale with wet feet. Do not stand on the scale wearing socks.
- The scales have a measuring range up to 180 kg (396 lb / 28 st 4 lb). Do not overload the scales – this can cause irreversible damage.
- Handle the scale with care. The scale must not be subjected to impacts or vibration. Do not fall or drop anything on the scales.
- Do not attempt to repair the device yourself in the event of a malfunction since this will invalidate the warranty. Repairs should only be carried out by authorised service centres.
- Children shall not play with the appliance.
- In the event of malfunctions, do not proceed with any unauthorised repair of the device. Otherwise, this will render any warranty claims void. Have repairs performed only by an authorised service centre.

WARNING – BATTERY SAFETY INFORMATION

- Do not take the batteries apart!
- If required, clean the battery and device contacts before insertion!
- Remove spent batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! In the event of contact with battery acid, immediately rinse the affected area with plenty of clean water and consult a doctor immediately!
- If a battery is swallowed, consult a doctor immediately!
- Insert the batteries correctly, paying attention to the polarity!
- Keep the battery compartment closed properly!
- If the device is not used for an extended period of time, remove the batteries!
- Keep the batteries out of the reach of children!
- Do not re-charge the batteries! There is a risk of explosion!
- Do not short-circuit the device! There is a risk of explosion!
- Do not throw into a fire! There is a risk of explosion!
- Store unused batteries in the packaging and not near metal objects to avoid a short circuit!

Voltage supply: 4.5V= (3x 1.5 V, AAA)

Disposal: This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Please remove the battery before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Batteries containing pollutant substances are marked as follows: Cd = contains cadmium, Hg = contains mercury, Pb = contains lead.

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