

IPX1

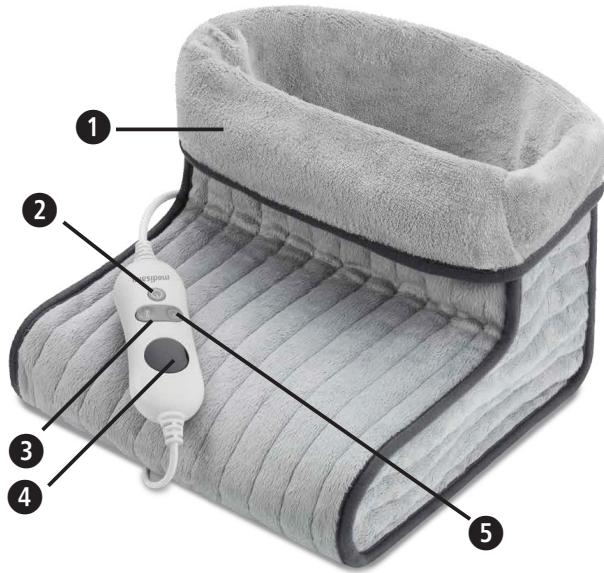Manual

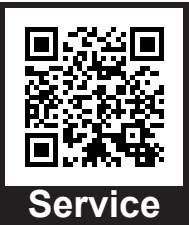Service

REF

EAN

medisana
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

60258

40 15588 60258 0

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Foot Warmer FWS

- 1 Cover (interior)
- 2 On/off key
- 3 Temperature level key
- 4 Display and indicator light
- 5 Operating time key

Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the packing is kept away from the reach of children! Risk of suffocation!
- Carefully check the foot warmer before each use for signs of wear and/or damage.
- Do not use if you notice wear, damage or signs of improper use on the foot warmer, the switch or the cables. Return it to the supplier.
- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- The appliance may not be switched on using a clock timer.
- The foot warmer may not be operated while folded or creased.
- No safety pins or other pointed or sharp objects may be attached to the foot warmer or stuck into it.
- The unit is designed for private home use only.
- This foot warmer is not intended for clinical or commercial use.
- Only use the foot warmer for the purpose described in this instruction manual.
- Do not use the device on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- The device must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians and unless the child has been well instructed how to use the control unit properly.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is not to be used for warming animals.
- Never leave the foot warmer unattended during operation.
- Do not fall asleep while the foot warmer is switched on.
- Do not wear shoes when using the foot warmer because this will lower your heat-sensitivity and will cause the foot warmer to become dirty.
- Prolonged use of the foot warmer at a high setting could lead to skin burns.
- Do not use the foot warmer in bed.
- Do not use the foot warmer when it is damp or in moist environments. The foot warmer may only be used again once it is completely dry.
- The controls and leads must not be exposed to any kind of moisture.
- Do not reach to retrieve a foot warmer that has fallen into water. Unplug the appliance immediately.
- During operation, the switch and lead may not be placed on or under the foot warmer or covered in any other way whatsoever.
- Do not place any heat sources (e.g. hot-water bottles or heating pads) on the foot warmer.
- Keep the mains lead away from hot surfaces.
- Do not place any large or heavy objects (e.g. cases, etc.) on the foot warmer regardless of whether the appliance is switched on or off.
- Never carry, pull or turn the foot warmer using the power cable and do not clamp the cable.
- When the appliance is not in use, disconnect the plug from the mains. In the event of a defect, do not repair the foot warmer yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- A damaged power supply cable cannot be replaced. If the power supply cable is damaged, the whole foot warmer must be disposed of.
- When storing the foot warmer, leave it to cool down first before folding. To prevent creases, do not place any objects on the foot warmer when storing.
- If you have any health concerns, consult your doctor before using the foot warmer.
- Do not use the foot warmer if your feet are swollen, inflamed or injured.
- If you experience any pain or discomfort while using the unit, stop using it immediately.

Voltage supply: 220 – 240V~ 50/60Hz max. 100W

Disposal: This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

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