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

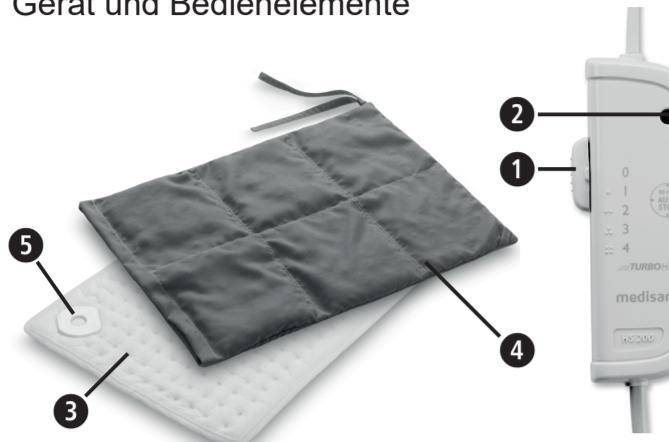
GB Instruction manual - Electric Cherry Stone Heat Pad HS 200

Device and controls

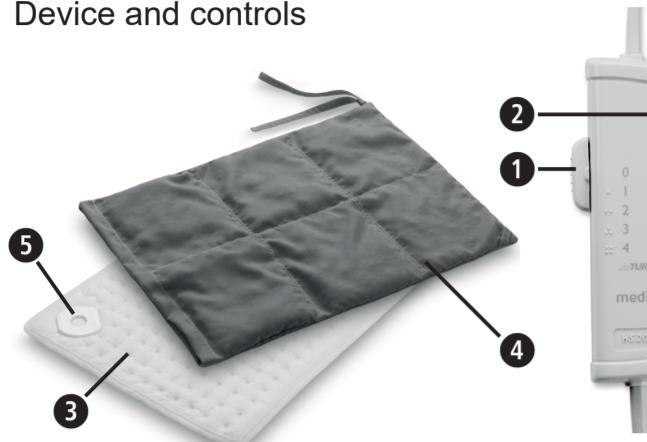
- ① Control unit with sliding switch
② Operating control display
③ Heat pad
④ Cherry stone cover
⑤ Plug connection at the back


Explanation of symbols

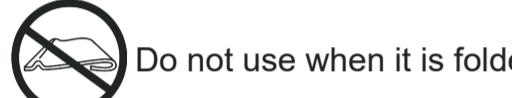Do not use when it is folded!

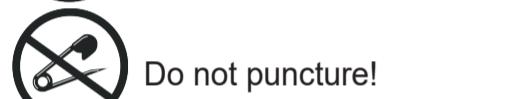
Do not puncture!

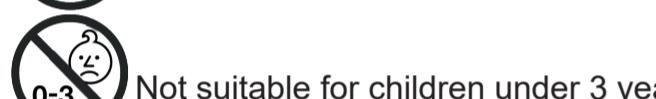
Not suitable for children under 3 years!


Only use indoors!

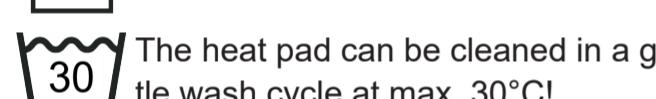
The heat pad can be cleaned in a gentle wash cycle at max. 30°C!

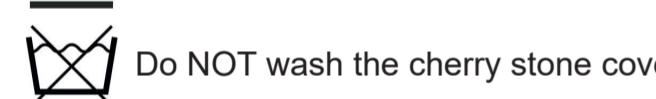
Do NOT wash the cherry stone cover!

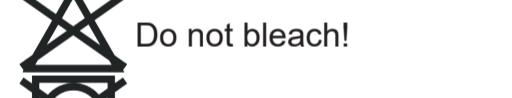Do not bleach!

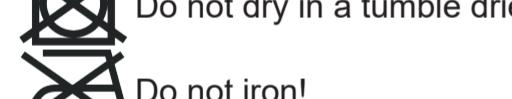
Do not dry in a tumble drier!

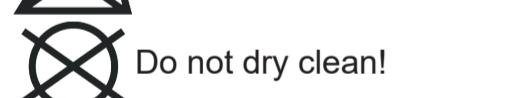
Do not iron!


Do not dry clean!

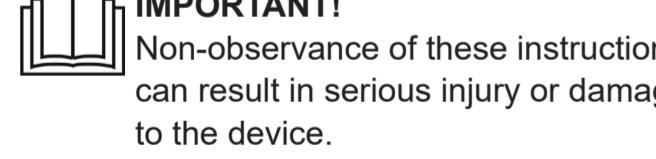**IMPORTANT!**


Non-observance of these instructions can result in serious injury or damage to the device.

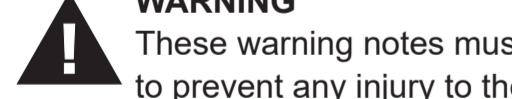**WARNING**


These warning notes must be observed to prevent any injury to the user.

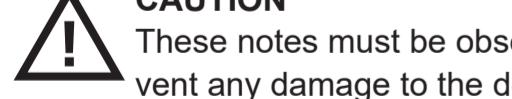**CAUTION**


These notes must be observed to prevent any damage to the device.

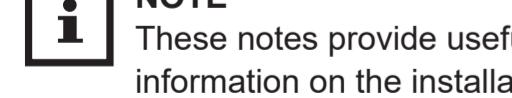**NOTE**


These notes provide useful additional information on the installation or operation.

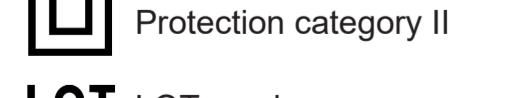Protection category II

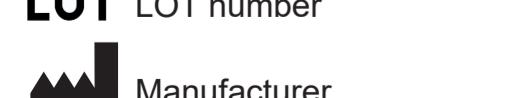Manufacturer

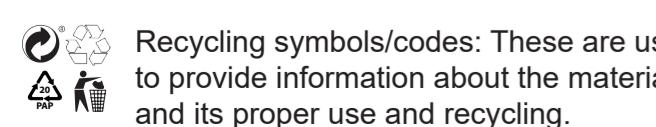GB Important instructions Retain for future use!

Power supply

- Before connecting the heat pad to the power supply, please check that it is switched off and that the mains voltage of your power supply corresponds to the voltage listed on the label.
- Keep mains cable and control unit away from heat, humidity and liquids. Never touch the mains plug when your hands are wet or when you are standing in water.
- Only operate the heating pad with the corresponding switching unit (HS 200).
- Do not attempt to touch a heat pad which has fallen into water. Pull the plug out of the socket immediately.
- To disconnect the heat pad, always pull the mains plug out of the mains socket. Never pull on the cable!
- Do not carry, pull or twist the heat pad using the mains cable.
- Make sure that the power socket and plug are easily accessible so you can unplug the mains quickly if necessary.

For special people


- Do not use the heat pad on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- The heat pad must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians and unless the child has been well instructed how to use the control unit properly.
- This appliance can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that electric cords shall be arranged in a safe way.
- Should you have any concerns about health issues, consult your doctor before using the heat pad.
- Do not use the heat pad on parts of the body that are swollen, injured or inflamed.
- Should you experience persistent pains in the muscles and joints, please inform your doctor. Persistent pains can be a symptom of serious illness.
- If you wear a pacemaker, please be aware that the electromagnetic fields emanating from this electric product may, under certain circumstances, interfere with the functions of your pacemaker. Therefore, ask your doctor and the manufacturer of the pacemaker for advice before using this product.
- If you experience any pain or discomfort while using the unit, stop using it immediately.
- Make sure your pets do not tear into your cherry stone cover; eating cherry pits can be poisonous to their little bodies.
- If the cover has come open and loose stones are leaking out, dispose of the product and the loose stones immediately and make sure they are out of reach of children and pets.
- The cherry stone cover is not suitable or allergic persons with appropriate allergies.

Operating the device

- Only use the device for its intended use as in the instruction manual. If misused, all guarantee rights shall become null and void.
- Carefully check the heat pad before each use, if you notice wear or damage.
- Do not use if you notice wear, damage or signs of improper use on the heat pad, the switch or the cables or if the device is not working.
- Do not fold or crease the heat pad during use.
- Neither pierce the heat pad with, nor attach safety pins or other spiky or sharp objects to the heat pad.

Risk of electric shock and fire! Never place the heat pad or controller in the microwave or on a heater.

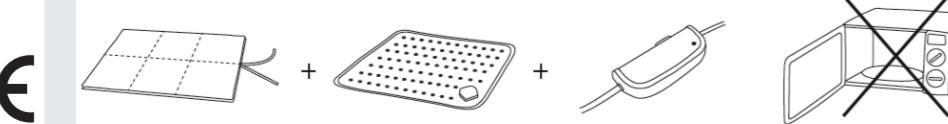Never use the heat pad unsupervised.

- The device is only for private use and is not to be used in hospitals.
- Do not fall asleep while the heat pad is switched on.
- Prolonged use of the heat pad at a high setting could lead to skin burns.
- Do not use the heat pad when wet and only use it in dry environments (not in the bathroom or similar environments).
- The control unit must not be placed on or under the heat pad during operation and it must not be covered.
- Do not sit on the heat pad, but place it like a cushion on the body part, as required.
- When using the cover as a heat pad, always remove the cover from heating pad.

Maintenance and cleaning

- You are only permitted to clean the unit. To avoid hazards, never repair the device yourself. Contact the customer service.
- Cleaning and user maintenance shall not be made by children without supervision.

If the power supply cable is damaged, it may only be replaced by **medisana**, an authorised dealer or by qualified personnel in order to avoid hazards.

- A damaged power supply cable must be replaced with a power cable of the same type.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- When storing the heat pad, allow it to cool down before laying out flat in a clean and dry place.
- Do not store the device with any objects covering or resting on top of it in order to avoid sharp wrinkles.

Scope of supply

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the it and contact your dealer or your service centre.

The following parts are included:

- 1 **medisana** Electric Cherry Stone Heat Pad HS 200 with removable switching unit HS 200
- 1 Instruction manual

If you notice any transport damage during unpacking, please contact your dealer without delay.

**WARNING**
Please ensure that the polythene packing is kept away from the reach of children!

Risk of suffocation!

Use


Do not sit on the heat pad, but place it like a cushion on the body part, as required. The heat pad reaches a pleasant temperature within a few minutes of switching on. The temperature can be adjusted as required.

In addition, the heat pad can also be operated in combination with the cherry pit cover. **Please note that it may take longer for the warmth to be felt on the outside of the cover.**

The cherry stone cover also offers the possibility of releasing the previously stored heat without a power source. To do this, warm up the cherry stone cover with the heat pad, then switch off the heat pad and separate the cover from the heat pad.

Separated from the heat pad, the *cherry stone cover* can also be warmed in the microwave (max. 700 watts, max. 1 minute) or on a heater. Spray around 6 ml of water evenly on the surface of the cover. Put 1/3 of a cup of water together with the cover into microwave. Always check the temperature carefully after heating to avoid the risk of burns.

CAUTION Do not use the hot cherry stone cushion together with the heat pad. Let the cherry stone cover cool down first (hand-warm) and only then use it as a cover for the heat pad.

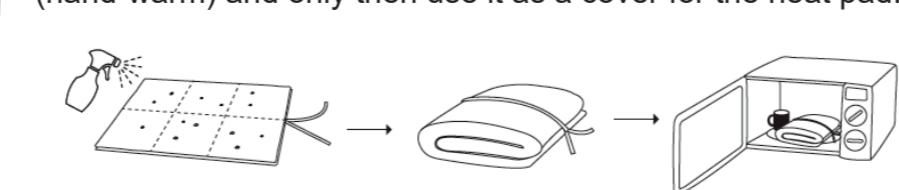Operation

Insert the mains plug into the socket and switch on the heat pad. To do this move the switch on the control unit ① to position 1, 2, 3 or 4 as required. The switch clicks into each position. When the appliance is switched on, the control display ② lights up and three short beeps are emitted. The higher the regulating level, the higher the temperature. To switch off the unit, move the switch completely back to position 0. If the operating control display goes out, the device is switched off. Move the selector switch to position 1 before falling asleep or any continuous use. The heat pad will switch off automatically after approx. 90 minutes of continuous use. Unplug the unit from the mains outlet once you have finished using it.

Switch positions

0 = the heat pad is switched off, the operating control display is off.

1 = lowest heat level, Level 1 is recommended for continuous operation.

2, 3 = middle heat level,

4 = highest heat level

The heat pad has an accurate electronic control system. Depending on the room temperature, the LED may change from orange to green colour and back within seconds. This indicates the fast re-heating and is normal.

Automatic cutout

The heat pad is equipped with an automatic cutout. This switches the unit off after 90 minutes of operation, independent of the heat level selected. Return the switch to position 0. However, to switch the heat pad back on, move the switch back to position 0 and then move it to the heat level required.

Cleaning and maintenance

Before cleaning the heat pad, unplug it from the mains outlet and let it cool down for at least ten minutes.

- The heat pad is equipped with a removable control unit. Unplug the control unit from the heat pad by pulling out the connection cable from the plug on the heat pad ⑤.
- The heat pad can be machine washed at 30° in a delicate programme or by hand gently. The best way to do this is to put the heat pad in a bath with lukewarm water and some mild detergent and then gently squeeze it out.

- Never use strong detergents or hard brushes.
- To dry, lay the heat pad flat on a moisture-absorbent towel or mat.

- Only use the heat pad again once it is completely dry.
- Reconnect the control unit to the heat pad in the correct manner.

- You should untangle the cable from time to time.
- Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.

- The cherry stone cover should be cleaned gently with a damp cloth. Do not wash it in the machine!

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Technical specifications

Name and model:	medisana Electric Cherry Stone Heat Pad HS 200
Power supply:	230V~, 50 Hz
Heat output:	100 watts
Automatic switch-off:	after approx. 90 minutes
Operating conditions:	Only use in dry rooms according to this user manual
Storage conditions:	Laid out flat and dry
Dimensions:	approx. 41 x 31 cm
Weight:	approx. 900 g
Mains lead length:	approx. 2.7 m
Art. number:	61168
EAN number:	40 15588 61168 1

In accordance with our policy of continual product improvement, we reserve the right for technical and design changes without notice.

The current version of this instruction manual can be found under www.medisana.com

Troubleshooting

This unit is equipped with the active protection system APS TECH. The unit detects different causes of errors and emits corresponding alarm sounds if necessary. In the case of an error, the power LED will light up in red.

3 beeps in 1-second interval and power LED lights up in red:
The connection cable of the control unit is not properly connected to the heat pad or another electrical defect is present. Please check the cable connection. Is the cable connected properly, please do not use the heat pad and contact the service center.

3 beeps in 2-second interval for one minute and power LED lights up in red:
The heat pad is broken, folded or fully covered. Please use the heat pad only according to this instruction manual and pay attention to all points in the safety information chapter.

Constant beeping and power LED lights up in red:
Important components are damaged or there is a short circuit. Please do not use the heat pad and contact the service center.

If there is no function at all after switching on the heat pad (power LED shows no light and no beeps can be heard), please check if the mains plug is correctly plugged in the mains outlet. If this is the case, the heat pad is defect. Please do not use the heat pad and contact the service center.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH

Carl-Schurz-Str. 2
41460 NEUSS
GERMANY

The service centre address is shown on the attached leaflet.

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