DE Massagesessel
GB Massage Chair
FR Fauteuil de massage
ES Sillón de masaje
NL Massagestoel

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GB Instruction Manual

1 Safety Information.....	6
2 Useful Information.....	8
3 Operating.....	8
4 Miscellaneous.....	9
5 Warranty.....	10

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
- GB**
- 1 Backrest
 - 2 Armrest
 - 3 Cable
 - 4 Control unit
 - 5 Pocket for control unit
 - 6 Seating area

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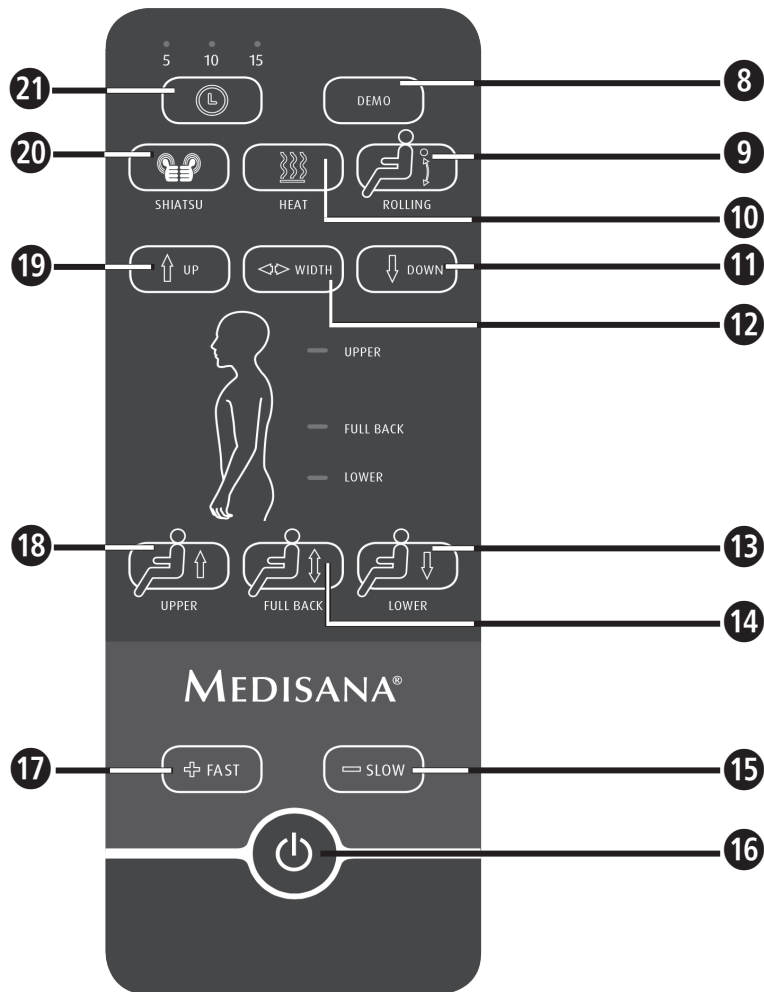- GB**
- 8 Demo button (performs all massage functions)
 - 9 Button for rolling massage
 - 10 HEAT button (switches the heating on/off)
 - 11 DOWN button (for positioning of the massage heads during shiatsu massage)
 - 12 WIDTH button (to change spacing of massage heads during rolling massage)
 - 13 LOWER button with LED (switches lower back massage on/off)
 - 14 FULL button with LED (switches whole back massage on/off)
 - 15 SLOW button (slower massage)
 - 16 On/Off button (switches the appliance on/off)
 - 17 FAST button (faster Massage)
 - 18 UPPER button with LED (switches upper back massage on/off)
 - 19 UP button (for positioning of the massage heads during shiatsu massage)
 - 20 Button for Shiatsu massage
 - 21 Timer button (for choosing massage time: 5, 10, 15 minutes)

Appareil et éléments de commande | Dispositivo y elementos de control | Toestel en bedieningselementen

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IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation.

Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

LOT LOT number

Manufacturer

**Protect from humidity!
Keep dry!**

Only to be used in enclosed spaces!

Protection category II

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Never operate the device with any power supply other than the original one.
- Do not operate in the vicinity of high-frequency electromagnetic transmitters.
- Keep the power cable and the unit away from heat, hot surfaces, moisture and liquids.
- Never touch the power plug or On/Off button with wet or damp hands or while standing in water.
- Unit components carrying voltage may not come into contact with liquid.
- Do not reach for the unit if it falls into water. Unplug the power cable immediately.
- Connect the unit so that the power plug is freely accessible.
- Always unplug the power plug from the power outlet immediately after use.
- Disconnect the unit from the mains by pulling on the plug. Never pull on the cable!
- Do not carry, pull or turn the unit by the power cable.

before using the unit

- Check the cable and the unit carefully before every use for damage. Do not put a defective unit into operation.
- Do not use the unit if damage to the unit or the cable is visible, if it does not function properly, if the cable or the unit have become wet.
- If the mains cable is damaged, it may only be replaced by **MEDISANA**, an authorized dealer or suitably qualified service personnel.

for using the unit

- Use the unit only for its intended purpose as according to the instruction manual. In case of misuse the warranty is no longer valid.
- Do not leave the unit unattended while it is connected to the mains.
- Avoid letting the unit come into contact with pointed or sharp objects.
- Never place the unit directly next to electric ovens or other heat sources.
- Do not stand on the unit.
- Only use this device in areas with sufficient space and where a good ventilation is guaranteed.
- Do not use the massage function for longer than 15 minutes. A too long massage can lead to overstimulation of muscles and to nerve irritations.
- Verify that seat supports and screw connections are in their correct place before each use.

**ATTENTION**

Do not move on the chair when the power cable is plugged in. This may damage the cable.

The unit is not intended for household or commercial purposes or for the medical sector. Consult your physician before using the chair if you have health concerns.

Do not use the massage chair in damp spaces (e.g. while bathing or showering).

- Do not expose the appliance to direct sunlight.
- To avoid the risk of electric shock, immediately pull out the plug from the power source:
 - if a malfunction occurs during operation
 - after each use
 - before each cleaning or maintenance.

for care and maintenance

- The chair is maintenance-free.
- In case of malfunction do not repair the unit yourself. This does not only annul the warranty, but can also cause serious danger (fire, electric shock, injury). Only let authorized service points carry out repairs.
- Cleaning and user maintenance shall not be made by children without supervision.

for special persons

- This unit is not intended for use by persons (including children below the age of 14) with restricted physical, sensory or mental abilities, lack of experience and/or lack of knowledge, unless they are supervised by someone responsible for their safety and have received instructions on how to use the unit.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit as an aid to or a replacement for medical applications. This may exacerbate chronic illness and symptoms.
- Do not treat any areas of the body, which show swellings, burns, inflammations, rashes, wounds or sensitive areas.
- You should refrain from using the device or consult your physician prior to use if:
 - you are pregnant,
 - you have a pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses or medical conditions: circulatory disorder, varicose veins, open wounds, contusions, skin cracks, phlebitis.
- The unit has a hot surface. Persons insensitive to heat must be careful when using the unit.
- Discontinue use and consult your physician if you feel pain or if the massage feels uncomfortable.

**WARNING**

Danger of tipping – Do not load the armrests or the back rest one-sided. Never let children use the device without close supervision!

Thank you Thank you for your confidence and congratulations!
In purchasing the **Massage Chair RC 420** you have acquired a quality product of **MEDISANA**. In order to achieve desired success and long-lasting pleasure with your **MEDISANA Massage Chair RC 420** we recommend that you carefully read the following tips for use and maintenance.

2.1 Items supplied and packa- ging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your supplier or your service centre.

The following parts are included:

- 1 **MEDISANA** Massage Chair **RC 420**
- 1 Power supply unit
- 1 Assembly manual
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your supplier without delay.

WARNING

**Please ensure that the polythene packing is kept away from the reach of children!
Risk of suffocation!**

ATTENTION

Please note the included assembly manual and be sure to tighten all screw connections.

3 Operating

The Massage Chair features three functions:

Fingertip massage (Shiatsu massage), rolling massage and warming. The warming function can only be used together with the shiatsu massage function in the back area. Do not use the unit for more than 15 minutes at a time. After approx. 15 minutes, the unit will switch off automatically. Always allow the unit to cool down fully before you use it again.

1. Connect the massage chair to the power supply unit and insert its plug into the mains socket. Ensure that the power supply is easily accessible.
2. Take a seat on the shiatsu massage seating area **6** and familiarise yourself with the functions.
3. Press the On/Off button **16** to switch on the device.

Switching off the device

Press the On/Off button **16**. The power LED begins to flash and the **massage heads travel back up to the park position**. When this position is reached, the LED extinguishes.

If the device is disconnected from power during operation, the massage heads will directly move to park position after reconnecting power.

3.1 Functions

Demo-Taste **8**

Press the DEMO button, the device starts automatically with introduction of different back setting with different functions. The LEDs show the current massage regions and functions. The device will stop once all functions worked.

Heat-Taste **10**

Press the heat button to activate the heat function; the LED on the control unit lights up red. The heat is radiated through the massage heads. Press button again to deactivate the heat function; the LED goes out.

Shiatsu massage 20**Lower button (lower back) 13**

The Shiatsu massage heads travel up and down the lower half of the back in a circular motion. The reversal of the direction takes place when the respective end positions are reached.

Upper button (upper back) 18

The Shiatsu massage heads travel up and down the upper half of the back in a circular motion. The reversal of the direction takes place when the respective end positions are reached.

Full button (full back) 14

The Shiatsu massage heads move from the "Park" position to the lower back, massaging the back in a circular motion whilst doing so. Once the undermost position has been reached, the direction of the massaging heads is reversed and the heads then travel up to the top of the backrest. Once the heads reach this point, the direction of rotation is reversed again and then travel down the back again.

UP/DOWN buttons 11/19 are used for individual positioning of the massage heads during shiatsu massage. Press and hold the button to move the massage heads to the position you like. The massage heads will stay there and rotate.

Rolling massage 9

Depending on your choice the massage heads move along the full back (Full button 14), the lower back (Lower button 13) or upper back (Upper button 18). In contrast to Shiatsu massage the massage heads will not rotate. Press WIDTH button 12 to change spacing of massage heads during rolling massage.

Press the buttons 20/9 again to deactivate the massage functions.

FAST /SLOW button 17/15

You can choose the speed of the massage with these two buttons: fast (+) and slow (-).

Timer button 21

With the timer button, you can determine the length of the massage. Press button 21 to switch the timer to 5 minutes (the LED "5" lights up). Pressing the button again switches the timer to 10 minutes (the LED "10" lights up). Pressing the button once more switches the timer to 15 minutes (the LED "15" lights up).

4.1 Cleaning and maintenance

- Before cleaning the device, make sure that it is turned off and the mains plug is removed from the socket. Allow the unit to cool down. Only clean the massage seat cover with a soft cloth slightly dampened with a mild detergent (e.g. soap suds).
- Under no circumstances should you use aggressive detergents, strong brushes, solvents, petrol, thinners or alcohol. These could damage the surface.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.

4.2 Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your local authority or your supplier for information about disposal.

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EAN number : 40 15588 88412 2

Loadable up to 150 kg

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under www.medisana.com

Warranty and repair terms

Please contact your supplier or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be corrected free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

 **MEDISANA AG**
Jagenbergstraße 19
41468 NEUSS
GERMANY

The service centre address is shown on the attached leaflet.

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