# medisana®

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### Explanation of symbols

#### IMPORTANT

Non-observance of these instructions can result in serious injury or damage to the device.

#### WARNING

These warning notes must be observed to prevent any injury to the user.

#### CAUTION

These notes must be observed to prevent any damage to the device.

#### NOTE

These notes give you useful additional information on the installation and operation.

Protection category II

LOT number

Manufacturer

Recycling symbols/codes: These are used to provide information about the material and its proper use and recycling.

## GB Safety instructions

**Read the instruction manual carefully before using this device, especially the safety instructions and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.**

- The device is intended for indoor private use only and should not be employed for commercial purposes.
- Only use the device for the according the instructions for use. The guarantee is invalidated if the device is misused.
- Avoid contact with sharp or pointed objects.
- Do not use this device if you are pregnant.
- Persons with a pacemaker, artificial joints or electronic implants should always consult a doctor before using the device.
- You must not use the device if you suffer from one or more of the following diseases or medical conditions: Circulatory problems, varicose veins, open wounds, diabetes, fever, heart problems, bruises, skin tears or other skin injuries, phlebitis, tumours, thrombosis, osteoporosis, spinal fractures or slipped discs.
- Do not use the appliance on the following parts of the body: Head, face, front of neck, chest, abdomen, bones or joints (e.g. elbows, knees, etc.), kidney and genital area.
- Excessive pressure when massaging with the device can lead to injuries!
- Special care is necessary when the device is used on, by or near children, sick and helpless persons.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- If you feel any pain from use of this device, or find the massage unpleasant, stop the treatment and consult your doctor.
- Keep the appliance away from heat, hot surfaces, moisture and liquids. Never touch the appliance with wet hands or when standing in water. Do not reach for a unit that has fallen into water.
- Do not cover the device when it is switched on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury.
- Never store or use the appliance near an electric stove or other heat sources.
- Do not sit or stand on the unit.
- Do not use the appliance together with other appliances.
- Do not use the appliance before going to bed. The massage has a stimulating effect and can delay falling asleep.
- Do not fall asleep while using the massager on yourself. Using it on the same spot for too long can cause discomfort and burns.
- Do not repair the device yourself in case of a malfunction. Repairs may only be carried out by an authorized dealer or an appropriately qualified person.
- Check the unit for damage before use. A defective unit must not be put into operation!
- Cleaning and user maintenance shall not be made by children without supervision.
- Switch off the device after each use.
- Do not use the device immediately after eating or drinking alcohol or when you are tired.
- Do not drop the appliance or subject it to strong shocks.
- Do not apply the appliance while operating vehicles or machinery.
- Make sure that the cable does not become a tripping hazard. It must not be kinked, pinched or twisted. If the cable is defective, it must not be used any longer.
- A damaged mains cable can be replaced with a mains cable of the same type.
- If you have long hair, it is advisable to fix the hair appropriately behind or on top of the head so that there is no risk of tangling when using the massager.

## Scope of delivery

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the appliance and contact your dealer or your service centre.

The following parts are included:

- 1 medisana Cellulite Massager AC 850
- 1 Adapter with connection cable
- 1 Instruction manual

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly. If you notice any transport damage when unpacking, please contact your dealer immediately.



#### WARNING

**Make sure that the packaging films do not get into the hands of children! There is a danger of suffocation!**

## Use

- Insert the connector plug of the adapter into the socket 3 on the product and then insert the plug into a socket.
- Slide your hand under the adjustable handle 1 - this is adjustable and can therefore be adapted to your hand.
- Gently press the unit with the massage rollers 4 against the part of the body you want to massage and first switch on the unit with the switch 2 to position 1 (low level).
- If you want the massage rollers 4 to rotate faster, move the switch 2 to position 2 (high level).
- Depending on how you feel, you can apply no pressure or light pressure. You can roll along different parts of the body, see „Application examples“.
- When massaging, do not concentrate on the same part of the body for too long (max. 1 minute), but change the massage area more frequently during the application.
- An application should not last longer than 15 minutes!
- Use the device regularly so that the positive effect will appear.

## Application examples

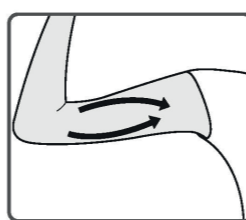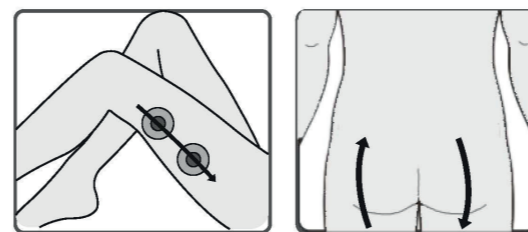Spot massage (hold the device on one spot)

Move the massager along

#### NOTES

- If you find direct application to the skin uncomfortable, you can also place a towel between the skin surface and the device during application.
- Application e.g. after bathing increases the circulation-promoting massage effect, as the body is already warm.

## Switching off the device

To switch off, slide the switch 2 back to the 0 position.

To switch off completely, disconnect the adapter from the mains socket.

## Cleaning and care

- For cleaning, make sure that the adapter with cable is not connected and let the unit cool down for at least 10 minutes after use.
- Never use aggressive cleaning agents, thinners, petrol, strong brushes or similar. The surface of the unit is sensitive.
- The unit should be cleaned with a soft, slightly damp cloth, possibly with a mild soap solution. To make cleaning easier, you can remove the massage rollers 4 and clean them separately. Do not mount the rollers again until they are completely dry.
- Clean the adaptor with cable with a clean, slightly moistened cloth.
- Never immerse the device or the adapter in liquids.
- Only use the unit again when it is completely dry and make sure that no moisture can get into the unit!
- It is best to store the unit in its original packaging and keep it in a clean and dry place out of direct sunlight.

## Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

## Technical specifications

Name and model: medisana Cellulite Massager AC 850

Power supply: Input: 100-240V~ 50/60Hz, 400mA  
Output: 12V=, 1200mA

Operating conditions: Only in dry rooms according to the instructions for use

Storage conditions: in a cool, dry place  
Dimensions: approx. 130 x 110 x 85 mm

Weight: approx. 220 g (without adapter)

Art. number: 88540

EAN number: 40 15588 88540 2

**In accordance with our policy of continual product improvement, we reserve the right for technical and design changes without notice.**

The current version of this instruction manual can be found under [www.medisana.com](http://www.medisana.com)

## Warranty and repair terms

Your statutory warranty rights are not restricted by our warranty described below. In the event of a warranty claim, please contact your specialist shop or the service centre directly. If you have to send in the unit, please state the defect and enclose a copy of the purchase receipt. The following warranty conditions apply:

1. The warranty period for medisana products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be corrected free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
  - a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.
  - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
  - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
  - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH  
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

You will find the service address on the separate enclosed sheet.

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