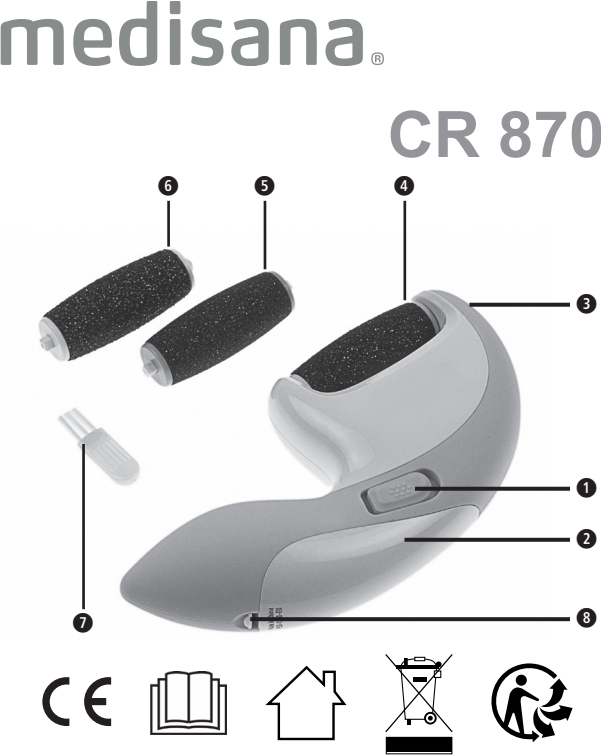
medisana
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

REF 88576

EAN 4015588 88576 1

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Electric callus remover CR 870

- 1 Power switch
- 2 Battery compartment (closure: ❸)
- 3 Release button
- 4 Fine roller (grey edge)
- 5 Coarse roller (blue edge)
- 6 Extra-coarse roller (white edge)
- 7 Cleaning brush

Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by **medisana**.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- You may only clean the device. In the event of malfunctions, do not proceed with any unauthorised repair of the device. Otherwise, this will render any warranty claims void. Ask your dealer and have repairs done only by **medisana** services.
- FOR HOUSEHOLD USE ONLY.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- If you are susceptible to allergic reactions or if you have circulatory problems, consult a physician before using this product.
- Do not use on irritated, sunburned or chapped skin, or on any area with dermatological problems.
- Do not apply over skin eruptions, warts, moles or varicose veins.
- Use should be easy and painless. If you experience nail or skin irritation, or hand discomfort when using this product, discontinue use. If a more extreme reaction occurs, consult your doctor.
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by Diabetics.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- Ensure that your hands are dry when operating the controls.
- This product should NEVER be used by any individual suffering from any ailments that would limit the user's capacity to operate the controls.
- Please ensure that all hair, clothing and jewellery are kept clear of moving parts of the product at all times.
- Before use, make sure that the roller is in good condition. Any malformation, distortion or damage may cause skin injury.
- Do not apply excess pressure or cause the motor to stop rotating. For best results, run the unit gently over the skin. This will ensure soft, safe and smooth removal of dead skin.
- Always turn the device off before cleaning it.
- Only clean the device with the supplied brush or a dry cloth. Do not use chemical agents including benzene, alcohol etc.
- Remove batteries if equipment is not going to be used for extended periods of time, and store batteries in a cool dry place away from direct sunlight and metal objects.
- If motor power seems reduced, replace the batteries.
- If roller surface becomes worn or less abrasive, fit a new roller.

BATTERY SAFETY INFORMATION

- Do not disassemble batteries!
- Clean the battery and device contacts if necessary before putting in the batteries!
- Remove discharged batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of this parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Keep the battery compartment well sealed!
- Remove the batteries from the device if it is not going to be used for an extended period!

- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Keep unused batteries in their packaging away from metal objects in order to prevent short circuiting!
- Do not throw used batteries into the household refuse; put them in a hazardous waste container or take them to a battery collection point, at the shop where they were purchased!

Inserting the batteries

- Open the battery compartment on the underside of the device.
- Insert batteries observing the polarity signs (+ and –) inside the battery compartment. NOTE: The batteries are stacked with one of the batteries behind the other two.
- Close the battery compartment.

How to use

- Ensure that your feet are clean and dry.
- Use the switch to select a power setting;
- For sensitive skin, or softer skin use the LO setting.
- For harder skin use the HI setting.
- Using light pressure place the surface of roller onto the heel or sole of your foot and gradually move the device around to target areas of hard skin. Do not use on a single area for more than 4 seconds at a time.
- Check the results throughout. Once you are happy with the results turn the switch to OFF.
- Rinse your feet with clean water & dry before applying your favourite foot cream.

Changing the rollers

To remove the roller.

- Ensure the switch is in the OFF position.
- With one hand press the release button on the side of the device.
- With the other hand lift out the roller. TIP: It's easier to lift the roller out from the end closest to the release button.
- To insert a new roller.**
- Loosely position one of the roller pins into the recess on the end furthest from the release button.
- Push the other end of the roller into place using the groves to guide it into position.
- The roller will secure into place with a click.

Voltage supply: 4.5 V ---
3 batteries 1,5 V --- , type AAA (micro), LR03

Disposal: This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Please remove the batteries before disposing of the device / unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your local authority or your supplier for information about disposal.

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